KXL® System

Cross-linking Solution For Your Keratoconus and Lasik Xtra® Patients







The World Leader in Corneal Remodeling™

KXL System

The new KXL System has a sleek, mobile design with a large touch screen which makes it easier to use.

Together with four preset protocols, a fifth user programmable treatment protocol and Avedro's family of high quality riboflavin formulations, you can choose a safe and predictable treatment for your patients via epi-off, epi-on, accelerated, pulsed or continuous treatment.

- Large 10.1" monitor for premium visibility and ease of use
- Live camera for improved treatment alignment and image acquisition
- On-board patient database for easy reference and planning
- Detailed patient treatment report for record keeping
- Four preset treatment protocols for speed and efficiency
- Fifth user programmable treatment protocol for convenience
- Intuitive software user interface








Treatment Presets For All Your Cross-linking Needs

TREATMENT PRESET	APPLICATION	FORMULATION
Epi-Off Accelerated Cross-linking	Cross-linking for keratoconus at high power and pulsed mode	Vibex Rapid
Epi-On Accelerated Cross-linking	Cross-linking for keratoconus at high power and pulsed mode without epithelial debridement	Trans-Epithelial Kit (Paracel)
Conventional Cross-linking	Dresden protocol cross-linking	Vibex Rapid
Lasik Xtra	Cross-linking after Lasik procedure after flap has been repositioned	Vibex Xtra
Programmable Treatment	Any user defined application	Any of the above as defined by user



Patient Database and Treatment Reports For More Efficient Patient Management

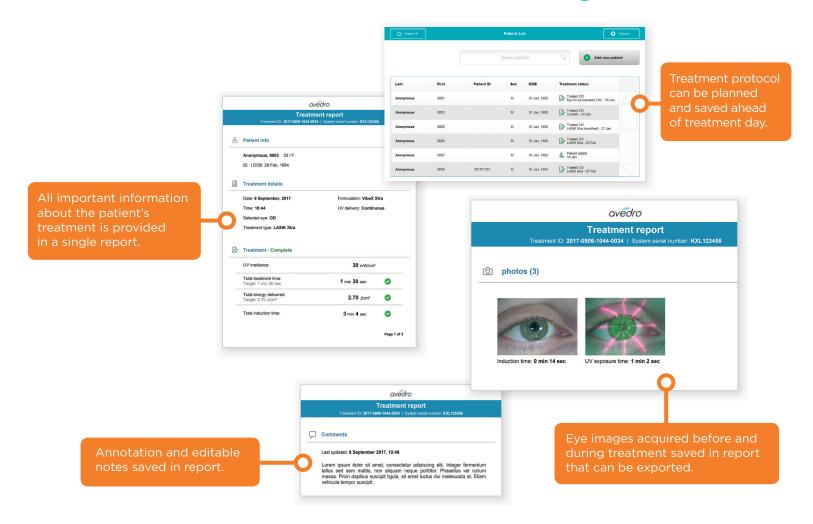

Cross-Linking Procedures Performed Across the Globe

