STELLAR^{M22} \

THE MULTI-APPLICATION PLATFORM YOUR PRACTICE

MUST HAVE

GROW YOUR PRACTICE.
MAXIMIZE YOUR REVENUE.



WWW.LUMENIS.COM

IMAGINE WHAT A SINGLE, MULTI-APPLICATION PLATFORM OFFERING 30 BESTSELLING, SKIN TREATMENTS AND HAIR REMOVAL COULD DO FOR YOUR PATIENTS AND YOUR PRACTICE...



THE PLATFORM THAT EXPANDS YOUR AESTHETIC OFFERINGS

With over 150 clinical studies, and Lumenis' 50+ years of innovation, this multi-application system offers over 30 cutting-edge skin and hair removal solutions for any gender, age and skin type.

Create Your Own Success Story and Elevate Your Practice with **CLINICALLY PROVEN RESULTS.**

REVENUE FROM DAY ONE

Real-world user data reveals it generates new revenue right from the start!

EXPAND YOUR BUSINESS

- ♦ Offer more treatments to more patients
- ♦ Become a one-stop solution for all your patients' aesthetic needs

GLOBAL TRUST

- → Recognized by practitioners worldwide for safety & efficacy
- ightharpoonup Leading the market with +30K systems sold in +100 countries

SCIENCE-BASED TREATMENT PRESETS

- → BASIC Mode with 40 preset protocols
- → NORMAL Mode with over 600 presets
- ADVANCED Mode for fully customized treatments

GROW YOUR PRACTICE OFFERING

4 TECHNOLOGIES IN 1 DEVICE



ONE VERSATILE APPLICATOR FOR 23 SKIN TREATMENTS AND HAIR REMOVAL

- 6th Generation IPL Device, invented by Lumenis in 1993.
- Optimal Pulse Technology (OPT®), controls pulse shape and duration.
- Multiple Sequential Pulsing (MSP™), enhances patient safety and comfort with cooling intervals between pulses
- 9 ExpertFilters™, easily interchangeable within seconds for versatile treatment options.
- 3 SapphireCool™ lightguides, ensure patient comfort with continuous contact cooling.



NEEDS JUST ONE PASS FOR EFFECTIVE SKIN RESURFACING, SAVING TIME AND PROTECTING PATIENTS' SKIN

- The only true fractional non-ablative technology.
- State-of-the-art CoolScan™ scanner with continuous contact cooling for increased patient comfort.
- Choose from over 600 combinations of area shape, size and density for optimal treatment, whether a single lesion or the whole face.
- Easy visibility of the treatment area through the SapphireCool™ Tip window.



SKIN TONING BY TREATING PIGMENTED LESIONS AND REMOVING DARK TATTOOS


- Homogenous 'Top Hat' beam profile for enhanced treatment safety and efficacy.
- Safe & Effective Treatments for All Skin Types Indicated for use in applications requiring selective photomechanical/photoacoustic effects.
- Single head for all Q-Switched applications.
- 7 spot sizes for customized treatments.



FOR PRECISE TREATMENT OF VASCULAR LESIONS AND LEG VEINS ON ALL SKIN TYPES

- Multiple Sequential Pulsing (MSP™) allows cooling between pulses, reducing risk of skin damage when treating Vascular Lesions. The series of sub-pulses enables deeper penetration and an overall longer pulse duration, which is crucial during treatment.
- 3 SapphireCool™ lightguides, changeable in seconds, for patient comfort through continuous contact cooling.

MORE TREATMENTS TO MORE PATIENTS

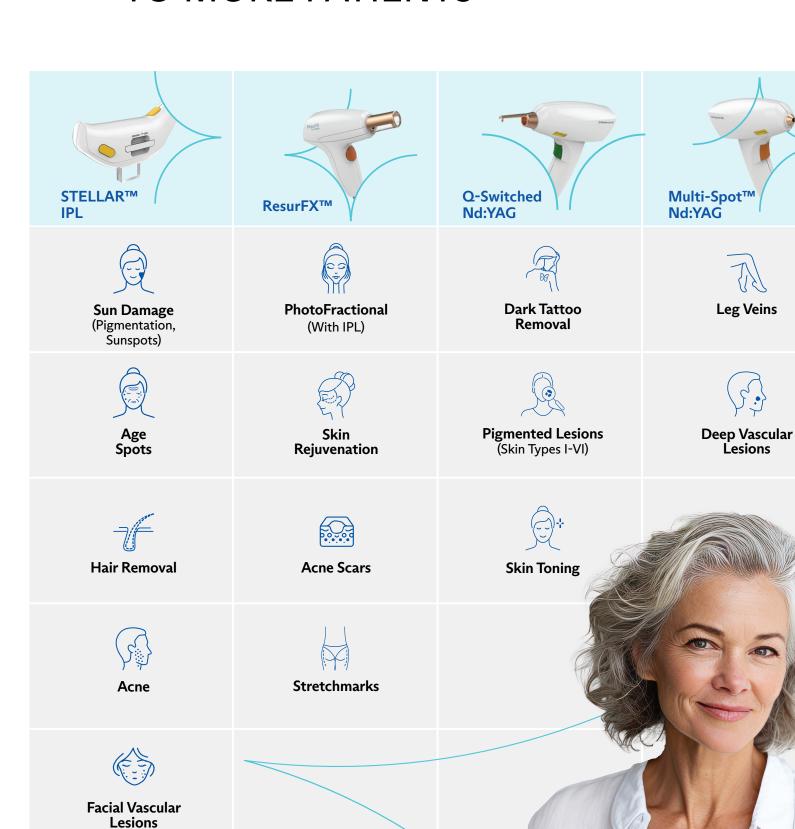

UNPRECEDENTED, LONG-LASTING RESULTS

NON-INVASIVE. TAILORED TO INDIVIDUAL NEEDS. MINIMUM SOCIAL DOWNTIME



ACNE TREATMENT WITH IPL Courtesy of Li Yuanhong, MD



SKIN RESURFACING WITH ResurFX™ Courtesy of Matteo Tretti Clementoni, MD



LEG VEINS TREATMENT WITH Nd:YAG Courtesy of Jose Miguel Montero Garcia, MD



Courtesy of Matteo Tretti Clementoni, MD



AFTER

ROSACEA TREATMENT WITH IPL Courtesy of Gilly Munavalli, MD





AFTER

Photofractional™ TREATMENT WITH IPL & ResurFX™ Courtesy of Matteo Tretti Clementoni, MD



TATTOO REMOVAL WITH Nd:YAG Q-Switched Courtesy of Michael Oeser, MD

When I think about purchasing a new device, I think about the ROI, and how quickly I'm going to be able to pay for the device... And I was happy to see that we paid for the device in just few months, which I think is unbelievable!

Melanie D. Palm, MD, MBA

STELLAR M22 IPL

9 EXPERT FILTERS | 3 LIGHTGUIDES **UNLIMITED POSSIBILITIES.**







Vascularity, IPL **Skin Treatments**



590nm Vascularity, IPL Skin Treatments, Pigmentation, Facial Veins



640nm Hair Removal, Lighter Skin



695nm Hair Removal, Darker Skin



755nm Hair Removal. Darker Skin

Pigmentation

Pigmentation,

Larger Facial Veins

615nm



Acne Vulgaris 400-600nm & 800-1200nm



Vascular Lesions 530-650nm & 900-1200nm



6 mm (Round)



8X15 mm



15X35 mm

CLINICAL INDICATIONS

- Poikiloderma of Civatte
- Angiomas Benign Epidermal Lesions
- Dyschromia Cutaneous Lesions
- Melasma Facial Telangiectasias
- Leg Telangiectasias
- Ephelides (Freckles)
- Erythema of Rosacea
- Hyperpigmentation
- Hemangiomas Leg Veins

Scars

- Hair Removal Port Wine Stains
 - Warts
 - Striae (Stretch Marks)
- Benign Cutaneous

Venous Malformations

Spider Angiomas

Vascular Lesions Truncal Telangiectasias

OPTIMAL PULSE TECHNOLOGY (OPT®)

Safe Energy Delivery:

OPT™ equalizes energy across the pulse, reducing the risk of burns and ensuring safer treatment.

• Enhance Treatment Comfort:

More controlled energy distribution leads to gentler, less painful procedures.


• Achieve Superior Results:

Consistent energy delivery improves efficacy, providing better, more predictable treatment outcomes.

IPL with OPT™

 Stellar M22 IPL with OPT. High-powered IPL without OPT.

 Limited capabilities of lower-powered systems

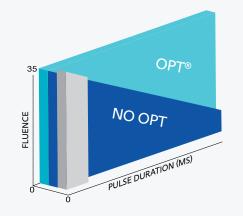

6 | STE*LL*AR^{M22}♦




- **♦ PRACTITIONER-FOCUSED DESIGN**
 - Clear visibility of treatment site using SapphireCool[™] Tip for safety and precision.


ResurFX[™] CoolScan[™] SHAPES:

6 Different Shapes, with Sizes and Densities changeable from the user interface, allow treatment precision and control.

No consumables, enhancing

practice profitability.

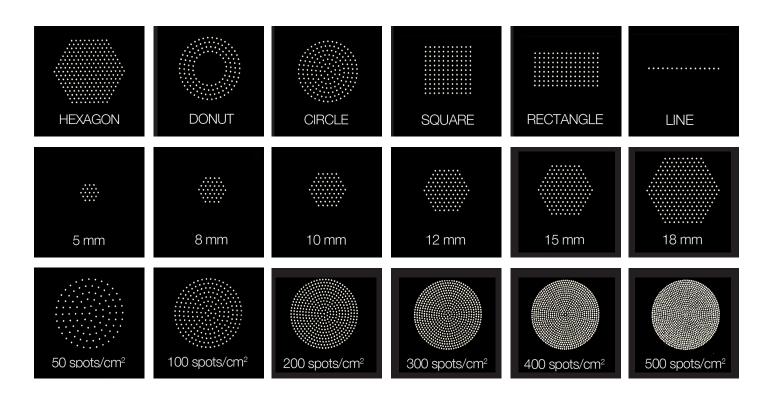

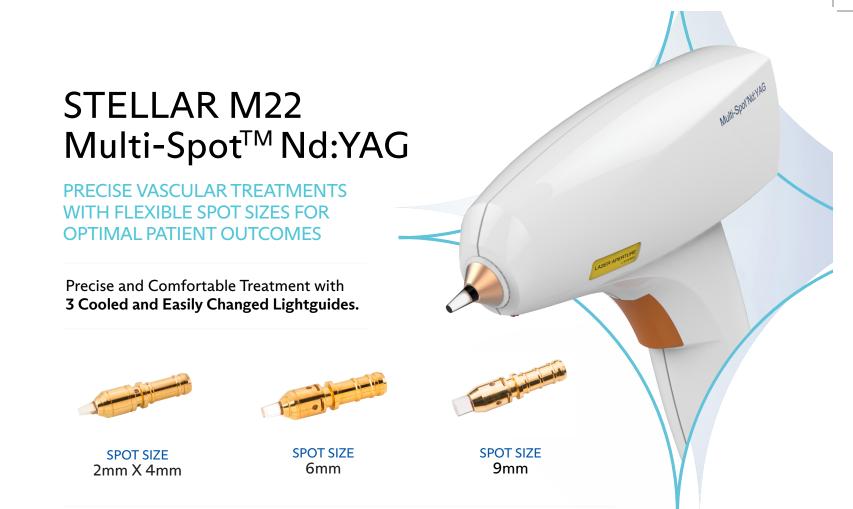

Maximized Treatment Safety with MSP™ Technology and preset parameters to ensure safe, consistent treatments every time.

MULTIPLE-SEQUENTIAL PULSING (MSP™)

• Deliver More Energy Safer:

MSPTM breaks the energy into multiple pulses, allowing surrounding tissue to cool between pulses. This ensures safer, depper, higher energy delivery without risking skin damage.

• Enhance Treatment Comfort:


By cooling the tissue between light pulses, MSPTM technology minimizes discomfort during treatments, creating a more pleasant experience for your clients.

• Provide Consistently Superior Results:


Percise energy control ensures uniform treatment, leading to more consistent and longer-lasting results.

→ Improves the efficacy in clearance of vascular lesions*

^{*} Pavan K.Nootheti, MD; Kimberly A. Pettit, MD; Gail Yosowitz, BS; Mitchel P. Goldman, MD. Clinical Improvement of Photodamaged Skin After a Single Intense Pulsed Light Treatment. The American Journal of Cosmetic Surgery Vol. 24, No. 1, 2007.

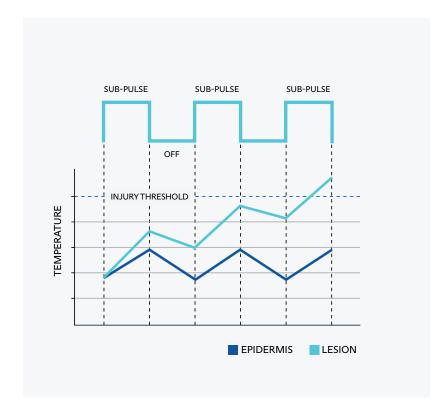


8 | STELLAR™22♦

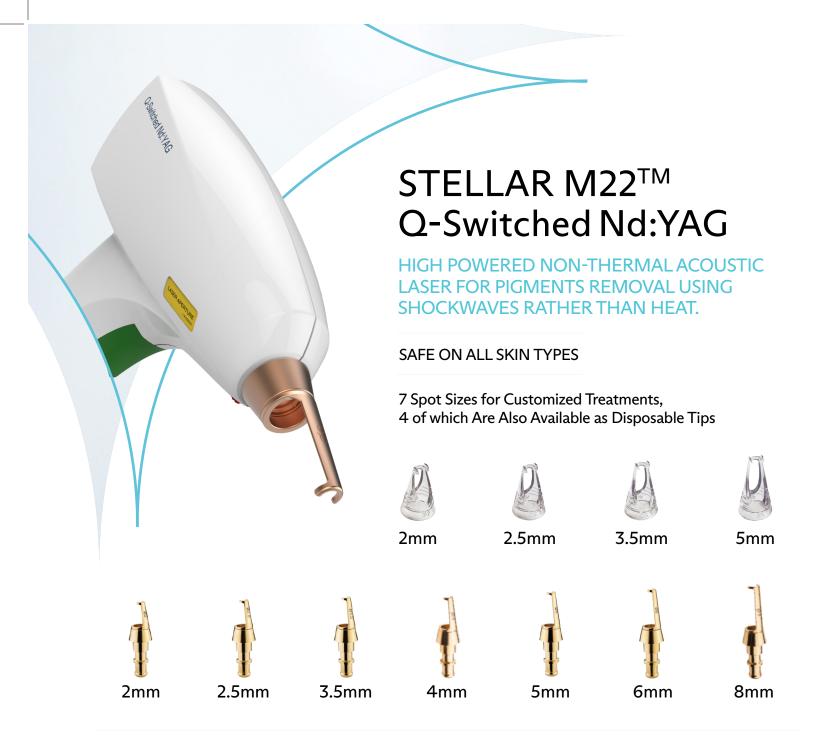


"TOP-HAT" BEAM PROFILE ADVANTAGES

- ♦ Enhanced treatment safety due to avoidance of bulk heating using the homogenizer.
- → Top-hat beam profile minimizes epidermal damage, tissue textural changes and the occurrence of scarring.
- ♦ A lack of a homogenizer might result in hot spots within the beam and increased risk of scarring.

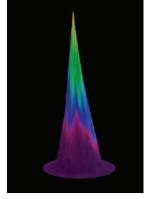

Beam profile side view without a homogenizer

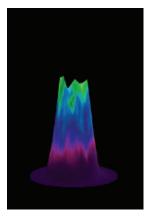

Beam profile side view with an homogenizer for a top-hat beam

YOUR JOURNEY WITH LUMENIS HAS JUST BEGUN

We want your business to grow with Lumenis. That's why we have put together a support package that begins with the first conversation and will stay with you as your business begins to see the benefits. We provide a support package which ensures that you have all the training you need and receive the most up to date clinical materials. Our direct to consumer marketing activities support the brand, keeping Stellar M22 the treatment of choice when it comes to laser and light skin treatment.



LUMENIS PRACTICE LOCATOR

Thousands of patients search for treatment providers every month. Help them find your clinic by registering your practice on the Lumenis Practice Locator. We run campaigns through influencers, social media and viral videos to promote this directory.



LUMENIS ONLINE

Consumers have many questions about laser hair removal treatments. Our dedicated patient-facing website will give them the re-assurance and confidence they need to move forward with the treatments.



TEAM TRAINING

We install the system and train your staff, so that you are able to assess skin types and conditions, creating the perfect treatment plan for each individual patient. Thanks to the clinical training you will be able to tailor treatments with confidence.



EXTENSIVE SERVICE NETWORK

We are committed to providing high quality, comprehensive service for medical equipment repair and maintenance. Lumenis service provides you with peace of mind, so you can focus on patient care, rather then on the equipment.



BRAND YOUR CLINIC

As a Lumenis client you will have access to the Lumenis Physician Portal. Here, you can find materials such as brochures and posters. As well as gain access to a wealth of clinical and product information.



YOUR MOST VALUABLE RESOURCE

Your dedicated representative is an email, text, or phone call away. We want you to feel part of the Lumenis family, so please feel free to contact us any time. We'll always do our best to get back to you asap.

SPECIFICATIONS

Stellar IPLTM

400-1200 nm
Acne (400-600 & 800-1200); Vascular (530-650 & 900-1200); 515 nm, 560 nm; 590 nm; 615 nm; 640 nm; 695 nm; 755 nm;
up to 56 J/cm2
4-20 ms
5-150 ms
Multiple Sequential Pulsing
Up to 1 Hz
35×15 mm²; 15×8 mm², 6mm
Continuous contact cooling

MultiSpot™ Nd:YAG

1064 nm
10-225 J/cm ²
2-20 ms
5-100 ms
Multiple Sequential Pulsing
Up to 1 Hz
6 mm; 2×4 mm ² ; 9 mm
Continuous contact cooling

ResurFX[™]

WAVELENGTH	1565 nm
PULSE ENERGY	10-70 mJ per micro-beam
PULSE CHARACTERISTIC	CoolScan™ Scanner
PATTERN SHAPES	Horizontal line, vertical line, square, horizontal rectangle, vertical rectangle, circle, donut, hexagon
BEAM DENSITY	50-500 micro-beams/cm²
REPETITION RATE	0.5-2 Hz
TIP	SapphireCool; Precision
SPOT SIZE	Up to 18 mm
COOLING	Continuous contact cooling

Q-Switched Nd:YAG

WAVELENGTH	1064 nm
FLUENCE	0.9-14 J/cm ²
PULSE DURATION	6-8 ns
PULSE CHARACTERISTIC	Top-Hat
REPETITION RATE	0.5-5.0 Hz
SPOT SIZE	2; 2.5; 3.5; 4; 5; 6 and 8 mm

SYSTEM SPECIFICATIONS

DIMENSIONS [without ResurFX]	51X56X52 cm (21X22X20.5 inches)	
DIMENSIONS [with ResurFX]	51X56X65 cm(21X22X25.6 inches)	
WEIGHT [without ResurFX]	45 kg / 99.2 lbs	
WEIGHT [with ResurFX]	60 kg / 132.3 lbs	
ELECTRICAL REQUIREMENTS	100 - 240 VAC ±10%, 12A max, 50/60 Hz, single phase, dedicated line	

WARNINGS AND RISKS

As with any laser or IPL, appropriate care must be taken to ensure safe and proper use.

The Stellar M22 operator's manual should be thoroughly reviewed and understood before operating this instrument. Stellar M22, including the Stellar IPL, Multi-Spot Nd:YAG, Q-Switched Nd:YAG and ResurFX, is contraindicated for patients with active infection, viral, fungal or bacterial diseases. It is also contraindicated for patients with skin disorders or inflammatory skin conditions. Be advised that IPL/Laser can cause epidermal injury. The risk increases with greater fluence/energy and skin pigmentation. The use of Stellar M22 could cause redness, swelling, change of pigmentation and scarring. See the system user manual for a complete list of contraindications and risk.

© 2025 Lumenis Be Ltd., All Rights Reserved. The Lumenis Group of Companies. Lumenis, its logo, Stellar M22, Stellar IPL, OPT, Multiple-Sequential Pulsing, ExpertFilters, SapphireCool lightguides, Photofractional, ResurFX and CoolScan are trademarks or registered trademarks of the Lumenis Group of Companies. Specifications are subject to change without notice Lumenis is ISO 13485:2016 certified.



Yokneam Industrial Park Hakidma Street,9. ISRAEL P.O.B 426 Yokneam 2069236

LUMENIS (Germany) GmbH

Heinrich-Hertz-Str 3 D-63303 Dreieich-Dreieichenhain Germany T +49 (0) 6103 8335 0

AMERICAS

San Jose, CA, USA

T +1 408 764 3000 +1 877 586 3647

F +1 408 764 3999

EMEA Dreieich, Germany

T +49 6103 8335 0 F +49 6103 8335 300

Formello (RM), Italy

T +39 06 90 75 230 F +39 06 90 75 269

Bucking hamshire, UK

T +44 20 8736 4110 F +44 20 8736 4119

Nice. France

T +33 4 8939 8341

ASIA / PACIFIC

Beijing, China

T + 86 1057 3766 77 T+82 2 6137 9225 F + 86 1057 3767 67 F+82 2 6137 9231

Seoul, Korea

T+61 1800 5863 64

F + 61 3 9863 7822

Gurgaon, India

Tokyo, Japan T +81 3 44318300 T+911244220795 F+911244854932 F +81 3 44318301

Kowloon, Hong Kong

T + 852 2174 2800 F + 852 2722 5151

Taipei, Taiwan

T+886 2 7719 9098 F: +886 2 7719 9005

WARNING





CE₀₄₈₃ Rx ONLY

© LUMENIS BE LTD. 2025

PB-00085660 Rev A

